




How to sew a surgical face mask

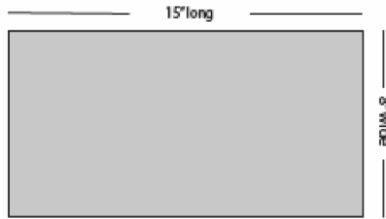


1

Cut one cotton rectangle, 8" wide by 15" long. Cut two 6.5" pieces of elastic.




2

Fold in half (right-sides together) and sew along the top edge, leave a 2" opening for the pocket.

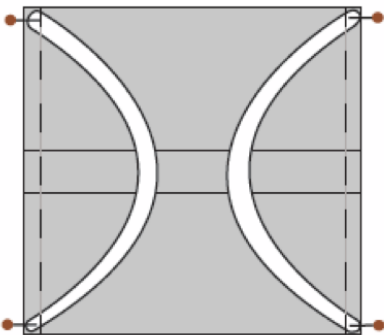

3

Sandwich the elastics between the layers of the mask. Pin in place at the corners, and sew along each side.

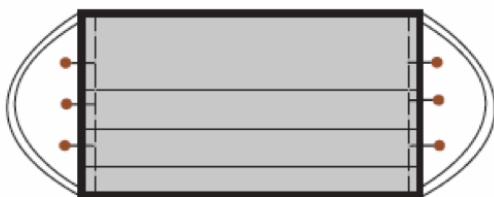

4

Turn the mask right side out. Fold 3 evenly spaced (about 1/2") pleats. Pin in place and sew to secure.

Any unauthorized copying, alteration, distribution, transmission, performance, display or other use of this material is prohibited